

Contents – Extractors

Nitinol retrieval devices

NGage® Nitinol Stone Extractor	5.2
NCompass® Nitinol Stone Extractor	5.3
NCircle® Nitinol Tipless Stone Extractor	5.4
NCircle® Delta Wire® Nitinol Tipless Stone Extractor	5.5
NCircle® Helical Nitinol Tipless Stone Extractor	5.6
NTrap® Stone Entrapment and Extraction Device	5.7
NForce® Nitinol Helical Stone Extractor	5.8
Perc NCircle® Nitinol Tipless Stone Extractor	5.9
Perc NGage® Nitinol Stone Extractor	5.10
Perc NCompass® Nitinol Stone Extractor	5.11

Stainless steel retrieval devices

Captura® Helical Stone Extractor	5.12
Atlas Wire™ Stone Extractor	5.13

Grasping forceps

Captura® Three-Prong Grasper	5.14
Grasping Forceps	5.15

Retrieval forceps

Alligator Tooth Retrieval Forceps	5.16
Rat Tooth Retrieval Forceps	5.17


Loops

NSnare™ Stent Retriever	5.18
-----------------------------------	------



NGage® Nitinol Stone Extractor

Used for urinary stone extraction. The patented design enables the physician to engage, reposition, release, or extract stones in the kidney or the ureter.

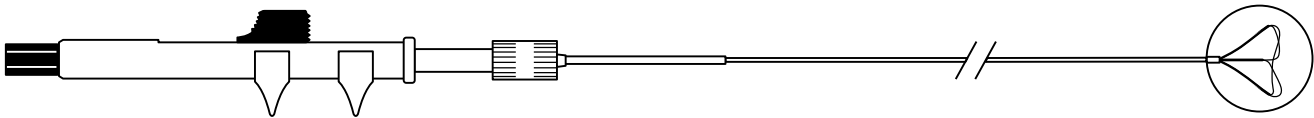

Order Number	Reference Part Number	Fr	Length cm	Expanded Diameter mm
G48294	NGE-017115	1.7	115	8
G48295	NGE-017115-MB	1.7	115	11
G48297	NGE-022115	2.2	115	8
G48296	NGE-022115-MB	2.2	115	11


Some products or part numbers may not be available in all markets. Contact your local Cook representative or Customer Support & Delivery for details.

Please see product risk information in the IFU at cookmedical.eu.



NCompass® Nitinol Stone Extractor

Used for urinary stone extraction. At full deployment, the multiwire geometry provides a traditional 3- or 4-wire basket. Partial closure redesigns the basket to provide the weave of a 12- or 16-wire basket. The duality of the NCompass extractor allows both routine capture and extraction of intact stones as well as the ability to capture multiple small fragments.

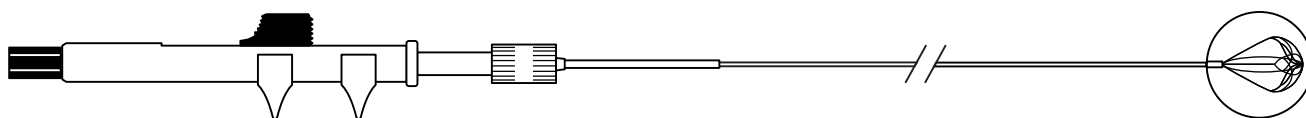

Order Number	Reference Part Number	Fr	Length cm	Basket Diameter cm	Basket Configuration
Tipless					
G46027	NCT4-017115	1.7	115	1.0	4 wire/16 wire
G32701	NCT4-024115	2.4	115	1.0	4 wire/16 wire
Modified Basket					
G32863	NCT4-024115-MB	2.4	115	1.5	4 wire/16 wire


Some products or part numbers may not be available in all markets. Contact your local Cook representative or Customer Support & Delivery for details.

Please see product risk information in the IFU at cookmedical.eu.



NCircle® Nitinol Tipless Stone Extractor

Used for stone manipulation and removal in the urinary tract. The nitinol wire allows the basket to retain its shape following extreme torsion. The tipless design allows the basket to open at the base of the calyx and encircle the stone. The tipless design also reduces the likelihood of mucosal damage and perforation.

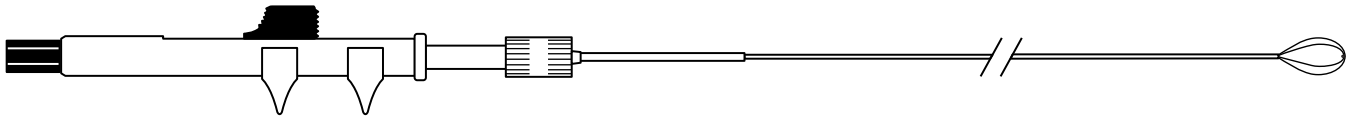

Order Number	Reference Part Number	Fr	Length cm	Basket Diameter cm
G46206	NTSE-015115	1.5	115	1.0
G18778	NTSE-022115-UDH	2.2	115	1.0
G17627	NTSE-030115-UDH	3.0	115	1.0
G17520	NTSE-045065-UDH	4.5	65	2.0
Modified Basket				
G18777	NTSE-022115-UDH-MB	2.2	115	2.0
G17872	NTSE-030115-UDH-MB	3.0	115	2.0


Some products or part numbers may not be available in all markets. Contact your local Cook representative or Customer Support & Delivery for details.

Please see product risk information in the IFU at cookmedical.eu.



NCircle® Delta Wire® Nitinol Tipless Stone Extractor

Used for stone manipulation and removal in the urinary tract. The Delta Wire increases the radial force of the basket. The tipless design allows the basket to open at the base of the calyx and encircle the stone. The tipless design also reduces the likelihood of mucosal damage and perforation.

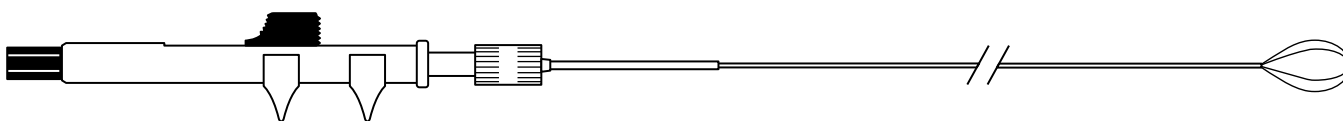

Order Number	Reference Part Number	Fr	Length cm	Basket Diameter cm
G19110	NTSED-024115-UDH	2.4	115	1.0
Modified Basket				
G19111	NTSED-024115-UDH-MB	2.4	115	2.0


Some products or part numbers may not be available in all markets. Contact your local Cook representative or Customer Support & Delivery for details.

Please see product risk information in the IFU at cookmedical.eu.



NCircle® Helical Nitinol Tipless Stone Extractor

Used for stone manipulation and removal in the urinary tract. The nitinol wire allows the basket to retain its shape following extreme torsion. The tipless design allows the basket to open at the base of the calyx and encircle the stone. The tipless design also reduces the likelihood of mucosal damage and perforation. The helical design aids in stone capture.

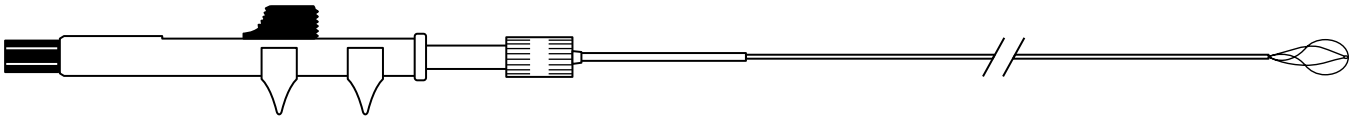

Order Number	Reference Part Number	Fr	Length cm	Basket Diameter cm
G18191	NTHSES-030115-UDH	3.0	115	1.0
G18193	NTHSE-045065-UDH	4.5	65	1.5


Some products or part numbers may not be available in all markets. Contact your local Cook representative or Customer Support & Delivery for details.

Please see product risk information in the IFU at cookmedical.eu.



NTrap® Stone Entrapment and Extraction Device

Used to entrap and extract calculi and other foreign bodies in the urinary tract, and to minimise stone migration during laser, ultrasonic, electrohydraulic, or pneumatic lithotripsy.

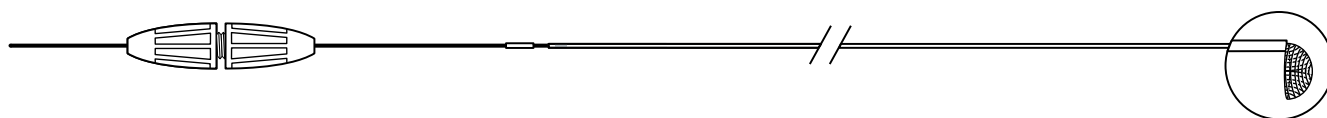

Order Number	Reference Part Number	Fr	Length cm	Basket Configuration
G32724	NTP-028145	2.8	145	7 mm "umbrella"


Some products or part numbers may not be available in all markets. Contact your local Cook representative or Customer Support & Delivery for details.

Please see product risk information in the IFU at cookmedical.eu.



NForce® Nitinol Helical Stone Extractor

Used for stone manipulation and removal in the urinary tract. The basket is constructed of nitinol Delta Wire®. This material allows the basket to retain its shape following torsion and increases the radial force of the basket. The helical design aids in stone capture.

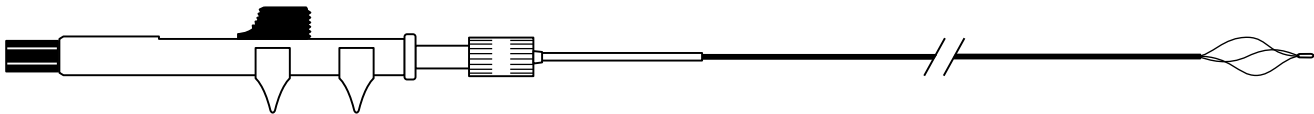

Order Number	Reference Part Number	Fr	Length cm	Basket Configuration
Occlusive Tip				
G19072	NFHSES-022115-3W-UDH	2.2	115	3 wire
G17940	NFHSES-032115-3W-UDH	3.2	115	3 wire


Some products or part numbers may not be available in all markets. Contact your local Cook representative or Customer Support & Delivery for details.

Please see product risk information in the IFU at cookmedical.eu.



Perc NCircle® Nitinol Tipless Stone Extractor

Used for the extraction of stones and foreign objects in the kidney and bladder. The tipless basket allows for direct positioning against the mucosal lining with minimal trauma, thus reducing bleeding and maintaining a clear working field. The nitinol wire construction allows the extractor to retain its shape following torsion. A stainless steel cannula provides strength and reliability comparable to traditional reusable surgical instrumentation.

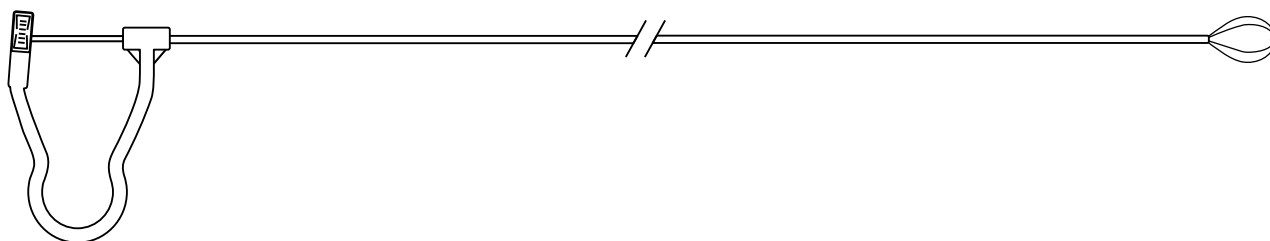

Order Number	Reference Part Number	Fr	Length cm	Basket Diameter cm
G32862	NTRSE-100038	10.0	38	2.0


Some products or part numbers may not be available in all markets. Contact your local Cook representative or Customer Support & Delivery for details.

Please see product risk information in the IFU at cookmedical.eu.



Perc NGage® Nitinol Stone Extractor

Used for the extraction of stones and foreign objects in the kidney and bladder. The Perc NGage holds like a basket and releases like a grasper. The patented design enables the physician to engage, reposition, release, or extract large calculi. The nitinol wire construction allows the extractor to retain its shape following torsion. A stainless steel cannula provides strength and reliability comparable to traditional reusable surgical instrumentation.

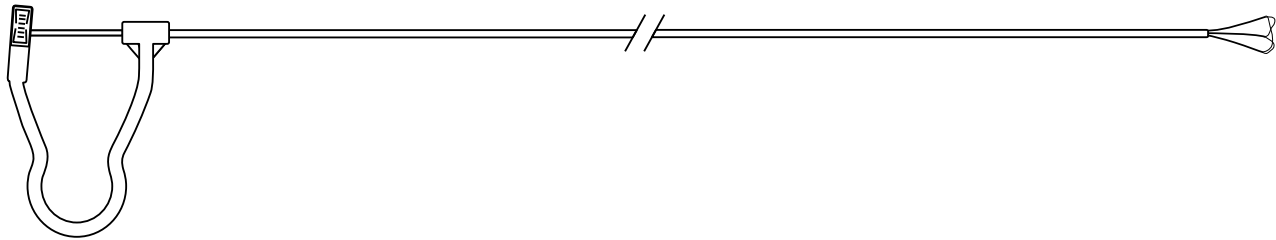

Order Number	Reference Part Number	Fr	Length cm	Basket Diameter cm
G56401	NGRE-100038	10.0	38	1.5


Some products or part numbers may not be available in all markets. Contact your local Cook representative or Customer Support & Delivery for details.

Please see product risk information in the IFU at cookmedical.eu.



Perc NCompass® Nitinol Stone Extractor

Used for the extraction of stones and foreign objects in the kidney and bladder. The multiwire geometry of the Perc NCompass creates a standard 4-wire design at complete deployment and an intricate 16-wire design during closure. At full deployment, the distal portion provides a netting effect to help prevent stone migration. The nitinol wire construction allows the extractor to retain its shape following torsion. The tipless basket minimises trauma to surrounding tissue. A stainless steel cannula provides strength and reliability comparable to traditional reusable surgical instrumentation.

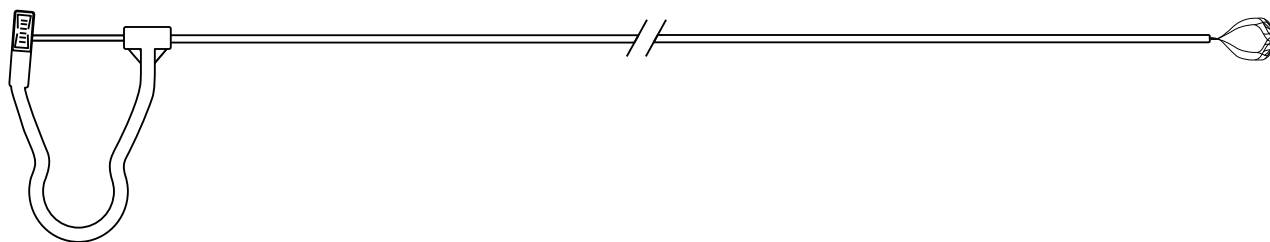

Order Number	Reference Part Number	Fr	Length cm	Basket Diameter cm
G34931	NCTR-100038	10.0	38	1.5


Some products or part numbers may not be available in all markets. Contact your local Cook representative or Customer Support & Delivery for details.

Please see product risk information in the IFU at cookmedical.eu.



Captura® Helical Stone Extractor

Used for stone manipulation and removal in the urinary tract. The helical design aids in stone capture.

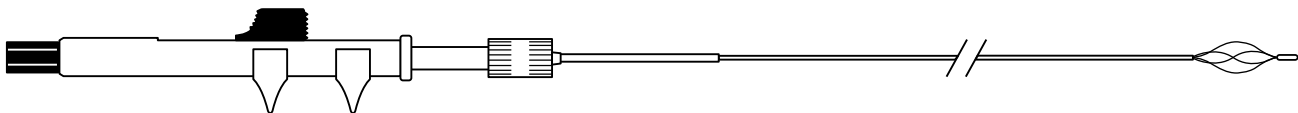

Order Number	Reference Part Number	Sheath Tip Fr	Sheath Fr	Length cm	Basket Configuration
Occlusive Tip					
G32757	CE-328115	2.8	2.5	115	3 wire
G32733	CE-332115	3.2	3.0	115	3 wire
G32756	CE-428115	2.8	2.5	115	4 wire
G32727	CE-432115	3.2	3.0	115	4 wire


Some products or part numbers may not be available in all markets. Contact your local Cook representative or Customer Support & Delivery for details.

Please see product risk information in the IFU at cookmedical.eu.



Atlas Wire™ Stone Extractor

Used for stone manipulation and removal in the urinary tract.

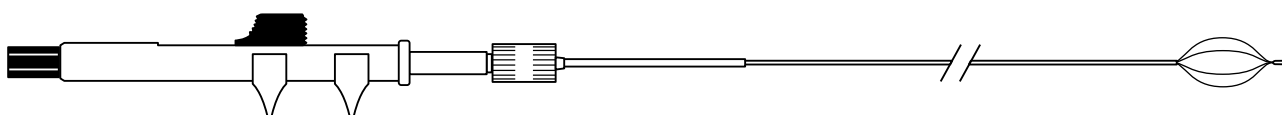

Order Number	Reference Part Number	Sheath Tip Fr	Sheath Fr	Length cm	Basket Configuration
Occlusive Tip					
G32734	AE-417115	1.7	1.7	115	4 wire
G32788	AE-428115	2.8	2.5	115	4 wire
G32824	AE-432115	3.2	3.0	115	4 wire


Some products or part numbers may not be available in all markets. Contact your local Cook representative or Customer Support & Delivery for details.

Please see product risk information in the IFU at cookmedical.eu.



Captura® Three-Prong Grasper


Used to facilitate the removal of stones and the retrieval of foreign objects under direct vision. Micro-throw grasping technology allows you to open the grasper wider in a shorter distance than standard graspers and to move the grasper in close approximation to the stone.

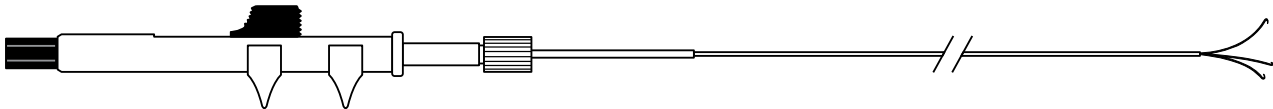

Order Number	Reference Part Number	Fr	Length cm	Max Grasping Diameter mm
G30361	112212-UDH-115-2.8	2.8	115	15
G30360	112212-UDH-115-3.2	3.2	115	15

Some products or part numbers may not be available in all markets. Contact your local Cook representative or Customer Support & Delivery for details.

Please see product risk information in the IFU at cookmedical.eu.

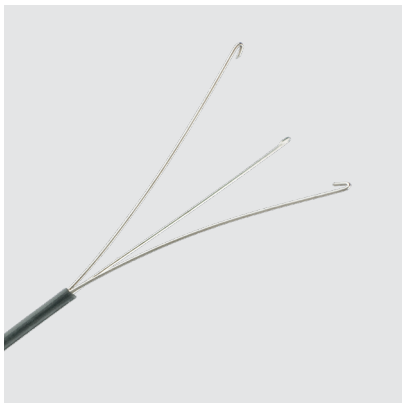

Grasping Forceps

Used for the extraction of stones and foreign objects in the urinary system.



Order Number	Reference Part Number	Fr	Length cm	Max Grasping Diameter mm	Comment
G16702	012212-UDH	2.5	115	10	with Unidex handle
G14596	012765	3.5	65	14	


Some products or part numbers may not be available in all markets. Contact your local Cook representative or Customer Support & Delivery for details.

Please see product risk information in the IFU at cookmedical.eu.



Alligator Tooth Retrieval Forceps

Used for the extraction of stones and foreign objects in the urinary system.

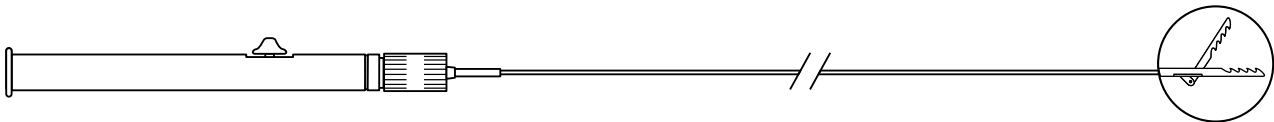

Order Number	Reference Part Number	Fr	Length cm
Single Jaw			
G14935	210321	3.3	65
G14922	210320	3.3	115


Some products or part numbers may not be available in all markets. Contact your local Cook representative or Customer Support & Delivery for details.

Please see product risk information in the IFU at cookmedical.eu.



Rat Tooth Retrieval Forceps

Used for the extraction of stones and foreign objects in the urinary system.

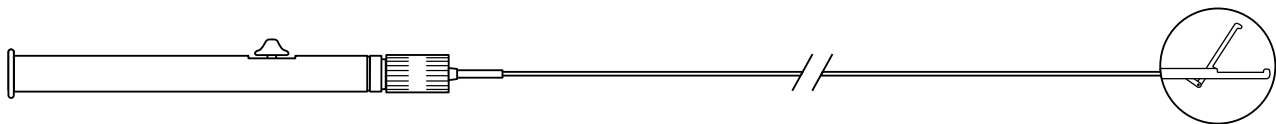

Order Number	Reference Part Number	Fr	Length cm
Single Jaw			
G14921	210310	3.3	115


Some products or part numbers may not be available in all markets. Contact your local Cook representative or Customer Support & Delivery for details.

Please see product risk information in the IFU at cookmedical.eu.



NSnare™ Stent Retriever

Used for the retrieval of stents and foreign objects in the urinary system. The flexible nitinol wire construction of the retrieval loop aids in stent capture and removal.

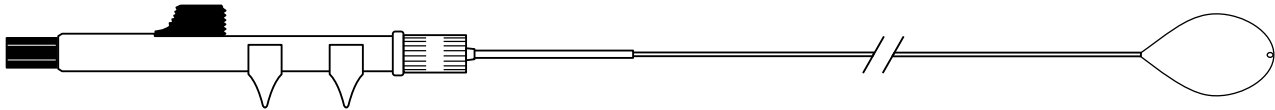

Order Number	Reference Part Number	Fr	Length cm
G46191	NSR-024115	2.4	115
G27233	NSR-045065-UDH	4.5	65

Some products or part numbers may not be available in all markets. Contact your local Cook representative or Customer Support & Delivery for details.

Please see product risk information in the IFU at cookmedical.eu.