




# ONE-STEP SUPRAPUBIC INTRODUCER



# ONE-STEP SUPRAPUBIC INTRODUCER

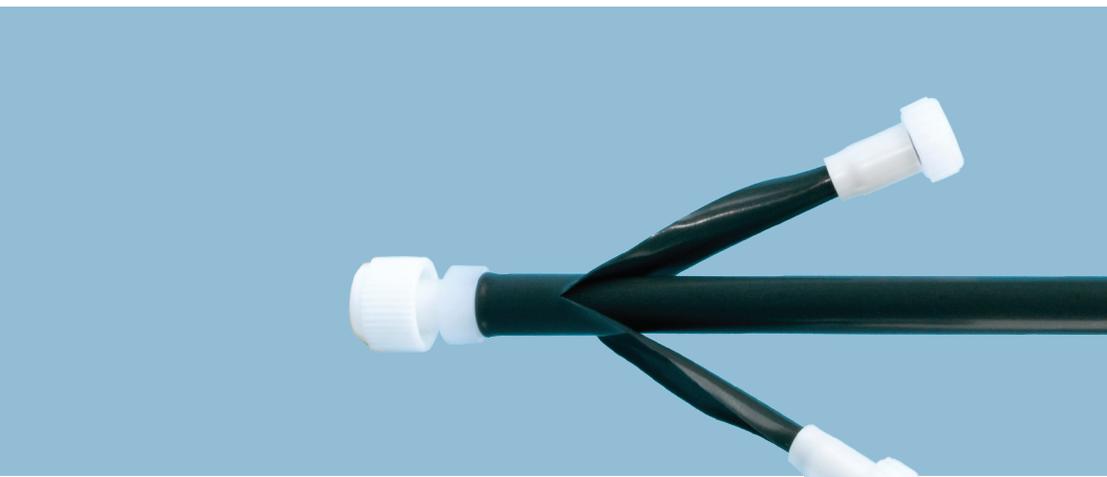

## The Convenience of One-Step

The One-Step Suprapubic Introducer provides access to place a drainage catheter suprapubically into the bladder. The removable trocar allows confirmation that the device is within the bladder by flow of urine. The dilator is removed to allow placement of a drainage catheter. The Peel-Away® sheath is then split into two parts, leaving the catheter in situ.

The One-Step Suprapubic Introducer provides the convenience of puncture, dilation and access in one simple step.

Set consists of trocar stylet, dilator and Peel-Away sheath. Supplied sterile in peel-open packages. Intended for one-time use.

Global Product Number	Order Number	Accepts Catheter* Fr	Sheath Length cm	Trocar Length cm
G25981	SANT-14-20-PLV	10	17	20
G26642	SANT-16-20-PLV	12	17	20
G26648	SANT-18-20-PLV	14	17	20
G26655	SANT-20-20-PLV	16	17	20
G26666	SANT-22-20-PLV	18	17	20

\*The outer diameter of catheters varies from manufacturer to manufacturer. Care should be taken to verify that the catheter will pass through the One-Step Peel-Away introducer prior to suprapubic insertion.



[www.cookmedical.com](http://www.cookmedical.com)

**COOK MEDICAL INCORPORATED**  
P.O. Box 4195, Bloomington, IN 47402-4195 USA  
Phone: 812.339.2235, Toll Free: 800.457.4500, Toll Free Fax: 800.554.8335

**COOK IRELAND LTD.**  
O'Halloran Road, National Technology Park,  
Limerick, IRELAND, Phone: +353 61 334440, Fax: +353 61 334441

**COOK (CANADA) INC.**  
111 Sandiford Drive, Stouffville, Ontario, L4A 7X5 CANADA  
Phone: 905.640.7110, Toll Free: 800.668.0300, Fax: 905.640.6804

**WILLIAM A. COOK AUSTRALIA PTY. LTD.**  
95 Brandl Street, Brisbane Technology Park, Eight Mile Plains  
Brisbane, QLD 4113 AUSTRALIA, Phone: +61 7 3841 1188

© COOK 2009 URO-BM-OSSIBR-EN-200906