



[www.cookmedical.com](http://www.cookmedical.com)

### Cook Europe Shared Service Centre

Country		Telephone	E-mail
Austria		(+43) 1-795-67-121	oeorders@cook.ie
Belgium	French	(+32) 27-00-16-33	beorders@cook.ie
Belgium	Flemish	(+32) 27-00-16-33	nlorders@cook.ie
Denmark		(+45) 38487607	daorders@cook.ie
Distributors		(+353) 61-239240	sscdistributors@cook.ie
France		(+33) 1-71-23-02-69	frorders@cook.ie
Germany		(+49) 6950-072804	deorders@cook.ie
Hungary		(+36) 17779199	huorders@cook.ie
Italy		(+39) 02-69-68-28-53	itorders@cook.ie
Netherlands		(+31) 202-013-367	nlorders@cook.ie
Norway		(+47) 23-16-2968	noorders@cook.ie
Poland		(+48) 223060159	plorders@cook.ie
Spain		(+34) 91-2702691	esorders@cook.ie
Switzerland	French	(+41) 44-8-009-609	frorders@cook.ie
Switzerland	Italian	(+41) 44-8-009-609	itorders@cook.ie
Switzerland	German	(+41) 44-8-009-609	deorders@cook.ie
Sweden		(+46) 85-876-94-68	seorders@cook.ie
UK		(+44) 20-7365-4183	ukorders@cook.ie

© COOK 2009

LIT0036-5

AORTIC  
INTERVENTION

CARDIOLOGY

CRITICAL  
CARE

ENDOSCOPY

INTERVENTIONAL  
RADIOLOGY

PERIPHERAL  
INTERVENTION

SURGERY

UROLOGY

WOMEN'S  
HEALTH



# Urology

PEDIATRIC FOCUS



# Contents

---

General Information.....	2
Injekt™ VUR Injection Needle .....	2
Pediatric Ureteral Catheter .....	3
Whistle Tip Ureteral Catheter .....	3
Angled Tip Ureteral Catheter .....	3
Flexor® DL .....	4
Flexor .....	4
Sof-Flex® Double Pigtail Stent Set.....	5
C-Flex™ Double Pigtail Ureteral Stent Set.....	5
Urethral Meatal Dilator .....	5
Tarkington Urethral Stent Set .....	6
Zaontz Urethral Stent .....	6
Koyle Diaper Stent .....	6
Mini-Perc™ Entry Set.....	7
Pediatric Nephrostomy Catheter/Stent Set .....	7
Salle Intraoperative Pyeloplasty Stent Set .....	7
Percutaneous Pigtail Nephrostomy Set .....	8
Silicone Foley Catheter .....	8
Pediatric Urinary Drain Bag Set .....	8

# General Information

---

## Pediatric Urology

---

Cook Medical is one of the world’s best-known and most respected names in the medical device industry. Closely allied with top medical research facilities, teaching hospitals and leading specialists around the world, Cook Medical continues to provide innovative, physician-driven, quality medical devices to healthcare providers who will accept nothing but the best for their patients.

Our pediatric urology products continue to provide cutting-edge technology to practitioners who demand better outcomes with fewer complications.

Cook Medical shares your high standards in pediatric patient care. We have provided this assurance for over 40 years and will continue to do so into the future.

We hope you’ll be pleased...

Please refer to our Focus Product Catalogue to view our range of Urology products.



### Injekt™ VUR Injection Needle

Used to inject legally marketed therapeutic/bulking agents into the genitourinary system, most specifically for VUR in the pediatric population and bladder neck injection in the general population.

Global Product Number	Order Number	Fr	Needle gage	Outer Sheath Length cm	Filiform Tip Length cm
G51158	IFI-383523	3.8	23	35	5



### Pediatric Ureteral Catheter

Used for drainage and retrograde pyelogram in pediatric patients. A stiffening stylet is included to aid in catheter placement.

Global Product Number	Order Number	Fr	Length cm
G15885	022230	3.0	70



### Whistle Tip Ureteral Catheter

Used for drainage and retrograde pyelogram. Sold in boxes of 10.

Global Product Number	Order Number	Fr	Length cm
<b>Right</b>			
G14028	022103	3.0	70
G14029	022104	4.0	70



### Angled Tip Ureteral Catheter

Used for directing a flexible-tipped wire guide in the ureter. The angled tip allows passage of the catheter past a ureteral stone in preparation for ESWL. Set includes wire guide.


Global Product Number	Order Number	Fr	Length cm	Wire Guide Diameter inch	Wire Guide Length cm	Tip Configuration
<b>20° Angled Tip</b>						
G14597	023104	4.0	70	.028	145	Tapered



**Flexor® DL**  
Dual Lumen Ureteral Access Sheath

Used to establish a conduit during endoscopic urological procedures facilitating the passage of endoscopes and other instruments into the urinary tract. The Flexor Dual Lumen (DL) Ureteral Access Sheath now provides ureteral dilation with two continuous working channels. The larger lumen provides protected introduction of endoscopes. The secondary lumen allows protected introduction of instruments 3 Fr or smaller, such as laser fibers, stone baskets or working wires. The additional lumen may be used for irrigation or the delivery of contrast media, freeing the working channel of the endoscope for additional uses.

Global Product Number	Order Number	Primary Lumen ID Fr	Secondary Lumen ID Fr	Length cm
G46072	FUS-095013-DL	9.5	3.0	13
G46029	FUS-120013-DL	12.0	3.0	13

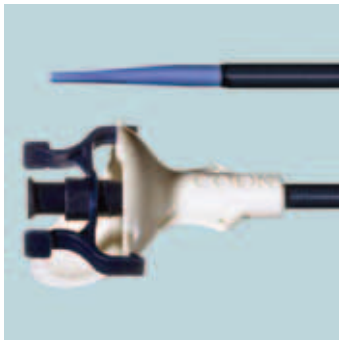

**Flexor®**  
Ureteral Access Sheath

Used to establish a conduit during endoscopic urological procedures facilitating the passage of endoscopes and other instruments into the urinary tract. The Flexor Ureteral Access Sheath provides ureteral dilation and a continuous working channel for the introduction of endoscopes and other instruments during ureteral access procedures. The Flexor sheath protects the ureter during repeated instrument exchanges, reducing the potential for trauma. The continuous channel also protects delicate instruments and the smaller flexible scopes from damage, reducing costly repairs.

Global Product Number	Order Number	ID Fr	Length cm
G32853	FUS-095013	9.5	13
G32852	FUS-120013	12.0	13



### Sof-Flex® Double Pigtail Stent Set

Used for temporary drainage from the ureteropelvic junction to the bladder in pediatric patients. Set includes stent, wire guide and stent positioner.

**CAUTION: Periodic evaluation is advised; stent must not remain indwelling more than 6 months.**

Global Product Number	Order Number	Fr	Length cm
<b>Set</b>			
G15606	039308	3.0	8
G15607	039310	3.0	10
G15608	039312	3.0	12
G15929	039314	3.0	14



### C-Flex™ Double Pigtail Ureteral Stent Set

Used for temporary internal drainage from the ureteropelvic junction to the bladder. Set includes stent, wire guide, stent positioner and catheter.

**CAUTION: Periodic evaluation is advised; stent must not remain indwelling more than 6 months.**

Global Product Number	Order Number	Fr	Length cm
<b>Set</b>			
G14794	036308	3.7	8
G14637	036310	3.7	10
G14834	036312	3.7	12
G14700	036314	3.7	14



### Urethral Meatal Dilator

Used for dilation of the penile urethral meatus.

Global Product Number	Order Number	Fr
G15161	073403	6.0–10.0



### Tarkington Urethral Stent Set

Used for stenting the urethra after hypospadias or epispadias repair and to provide postoperative drainage of the bladder. The stent is designed to be stabilized in the urethra without use of sutures at the meatus. Set includes stent and wire guide.


Global Product Number	Order Number	Fr	Length cm
G15614	032408	8.0	10



### Zaontz Urethral Stent

Used for stenting the urethra after hypospadias or epispadias repair and to provide postoperative drainage of the bladder.

Global Product Number	Order Number	Fr	Length cm
G16758	032306	6.0	12
G15049	032308	8.0	12
G15050	032310	10.0	12

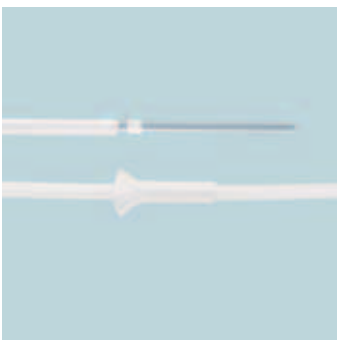

### Koyle Diaper Stent

Used for stenting the urethra after hypospadias or epispadias repair and to provide postoperative drainage of the bladder. The externalized portion of the stent is of sufficient length to allow drainage into a second or outside diaper, keeping the inside diaper dry during the healing process.

Global Product Number	Order Number	Fr	Length cm	Comment
<b>Stent only</b>				
G18827	KDSCF-060050-FL26-IM	6.0	50	26 cm from tip to flare with ink marks
<b>Set</b>				
G17666	KDSCFS-060050	6.0	50	
G17668	KDSCFS-080050	8.0	50	





### Mini-Perc™ Entry Set

#### Docimo Set

Used to gain percutaneous access into the renal pelvis to treat stone disease. The dual lumen catheter and the Peel-Away® introducer are sized for use in pediatric patients. Set includes needle, wire guides, sheath, introducer and catheter.

Global Product Number	Order Number
G17895	080000-DOCIMO



### Pediatric Nephrostomy Catheter/Stent Set

Used to prevent postoperative stricturing and to allow renal drainage after an endopyelotomy. **CAUTION: Periodic evaluation is advised; stent must not remain indwelling more than 4 months.**

Global Product Number	Order Number	Fr	Length cm	Stent Segment Fr	Stent Segment Length cm	Accepts Wire Guide Diameter inch/mm
G15034	080208-S6	8.2	15	4.7	5	.038/.97
G15035	080208-S7	8.2	15	4.7	7	.038/.97



### Salle Intraoperative Pyeloplasty Stent Set

Used to establish and maintain nephrostomy drainage following pyeloplasty. Set includes wire guide, stylet, catheter, stent, Tuohy-Borst adapter, retention disc and connecting tube. **CAUTION: Periodic evaluation is advised; stent must not remain indwelling more than 4 weeks.**

Global Product Number	Order Number	Fr	Usable Length cm	Wire Guide Diameter inch/mm	Wire Guide Length cm	Comment
G32773	SIPSF-040018-59	4.0	12-18	.025/.64	145	
G18168	SIPSF-040018-56-6	4.0	12-18	.025/.64	145	6 cm stylet
G32774	SIPSF-050018-59	4.7	12-18	.038/.97	145	



### Percutaneous Pigtail Nephrostomy Set

Used for percutaneous placement of a pigtail catheter in the renal pelvis for nephrostomy drainage. EchoTip® needle design provides enhanced visualization of needle tip when used with ultrasonic imaging equipment. Set includes two needles, wire guide, catheter, 7 dilators, connecting tube and retention disc.

**CAUTION: Periodic evaluation is advised; stent must not remain indwelling more than 4 months.**

Global Product Number	Order Number	Fr	Length cm
G14161	080006	6.0	15



### Silicone Foley Catheter

Used to provide drainage of urine from the bladder. All-silicone construction avoids latex sensitivity concerns. Sold in boxes of 12.

Global Product Number	Order Number	Fr	Length cm	Inflation Volume mL
<b>Open End</b>				
G17203	028506-OE	6.0	28	1.5
G17200	028508-OE	8.0	34	1.5
<b>Set with .021 inch Wire Guide (sold individually)</b>				
G17214	028506-OE-21	6.0	28	1.5



### Pediatric Urinary Drain Bag Set

Used to provide pediatric patients with an appropriately sized disposable system for urine drainage and collection. Supplied nonsterile. Intended for single-patient use only.

Global Product Number	Order Number	Volume mL
G17939	UDPS-200	200