




HiWire™ Nitinol Core Wire Guide

The HiWire™ is used to gain ureteral access in routine or demanding cases requiring precise control. Once in position, the HiWire may be used to facilitate instrument exchange or placement. Nitinol core allows for maximum deflection without kinking. Tapered tip designed to optimize flexibility. Available in standard and stiff shaft for use in navigating difficult anatomy, straightening of tortuous ureters, and facilitating instrument exchange or placement. 1:1 torque control allows for precise navigation throughout the most challenging anatomy. Hydrophilic coating for effortless advancement. Sold in boxes of five (5).

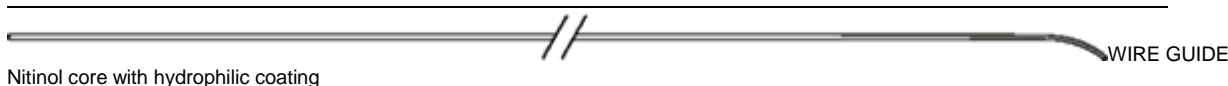

Global Product Number	Order Number	Diameter	Length	Tip Configuration	Tip Length
Standard Body Wire Guides					
NS	HW-025150	.025 inch	150 cm	Straight	3cm
NS	HWA-025150	.025 inch	150 cm	Angled	3cm
Stiff Body Wire Guides					
S	HWS-025150	.025 inch	150 cm	Straight	3 cm

U.S. Patent Number: 5,989,266

Kód VZP: 0092057

Mgr. Michaela Gray
Mob: 602 561 576